



THE QUEEN'S AWARDS  
FOR ENTERPRISE  
INTERNATIONAL TRADE  
2022

# PROTEIN RESIDUE TESTING

Sychem, in UK exclusive partnership with **Terragene**, provide a range of Protein Residue testing systems to a variety of industries.

MiniPro  
PRO1 MICRO  
Partnership  
PRO1 VT  
Bionova Cloud



Contact us today!

**0845 644 6824**



**hello@sychem.co.uk**



**sychem.co.uk**

Sychem and Terragene have worked in partnership since 2017.



## MiniPro

### Auto-reader Incubator


An Auto-Reader incubator for Hygiene Monitoring systems. The MiniPro has 3 positions for the incubation and transmittance readout of protein pens.

- 3 positions (2.5 cm depth, 1.3 cm diameter) for incubation
- Transmittance readout of hygiene indicators
- Results within 4 minutes
- Orifice for external thermometer
- Built-in top cover
- USB connection for recording results on PC using a Readout and Traceability software
- It allows to incubate and read 3 indicators simultaneously
- Thermal Printer results
- Sound alarm for indication of events.
- Firmware update by the user
- Dimensions: 102 mm high, 156 mm wide and 146 mm deep

## PRO 1 MICRO

High sensitivity protein residue testing pens. Fully quantitative and traceable results when used with the MiniPro Incubator.

- For the incubation and readout of Chemdye® PRO 1 MICRO
- Detection limit: 50 µg

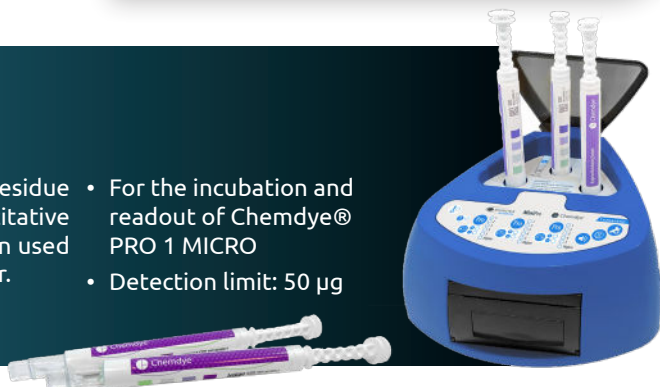



## UK Exclusive partnership with Terragene



Sychem are proud to be the UK exclusive partner for the Infection Prevention giant Terragene. Together we offer an extensive portfolio of products designed to fight against infection and tackle contamination.

Terragene offer a diverse range of Protein Residue Testing devices, including:

- MiniPro
- PRO 1 MICRO
- PRO1 VT



### Pro1 VT

A positive and negative validation kit to validate the efficacy of Chemdye® PRO1 MICRO Hygiene Monitoring System when used with Bionova® MiniPro.

The Negative Control consists of a vial with a protein-free solution (the green cap). The Positive Control consists of a vial with a defined amount of standard protein (the blue cap).



### Bionova Software

Terragene's technologically advanced system has the ability to quantify your results and keep all of your data in one place!

Protein Residue Testing made easy.



Contact us today!

**0845 644 6824**